




(86) Date de dépôt PCT/PCT Filing Date: 2004/09/22
 (87) Date publication PCT/PCT Publication Date: 2005/04/07
 (45) Date de délivrance/Issue Date: 2012/07/24
 (85) Entrée phase nationale/National Entry: 2005/06/14
 (86) N° demande PCT/PCT Application No.: US 2004/031202
 (87) N° publication PCT/PCT Publication No.: 2005/031356
 (30) Priorité/Priority: 2003/09/23 (US60/505,392)

(51) Cl.Int./Int.Cl. *G01N 33/558* (2006.01)
 (72) Inventeurs/Inventors:
 SMITH, HENRY J., US;
 ZAIDI, ASAD, US
 (73) Propriétaire/Owner:
 EPINEX DIAGNOSTICS, INC., US
 (74) Agent: SMART & BIGGAR

(54) Titre : TEST RAPIDE POUR ALBUMINE GLYCATEE
 (54) Title: RAPID TEST FOR GLYCATED ALBUMIN

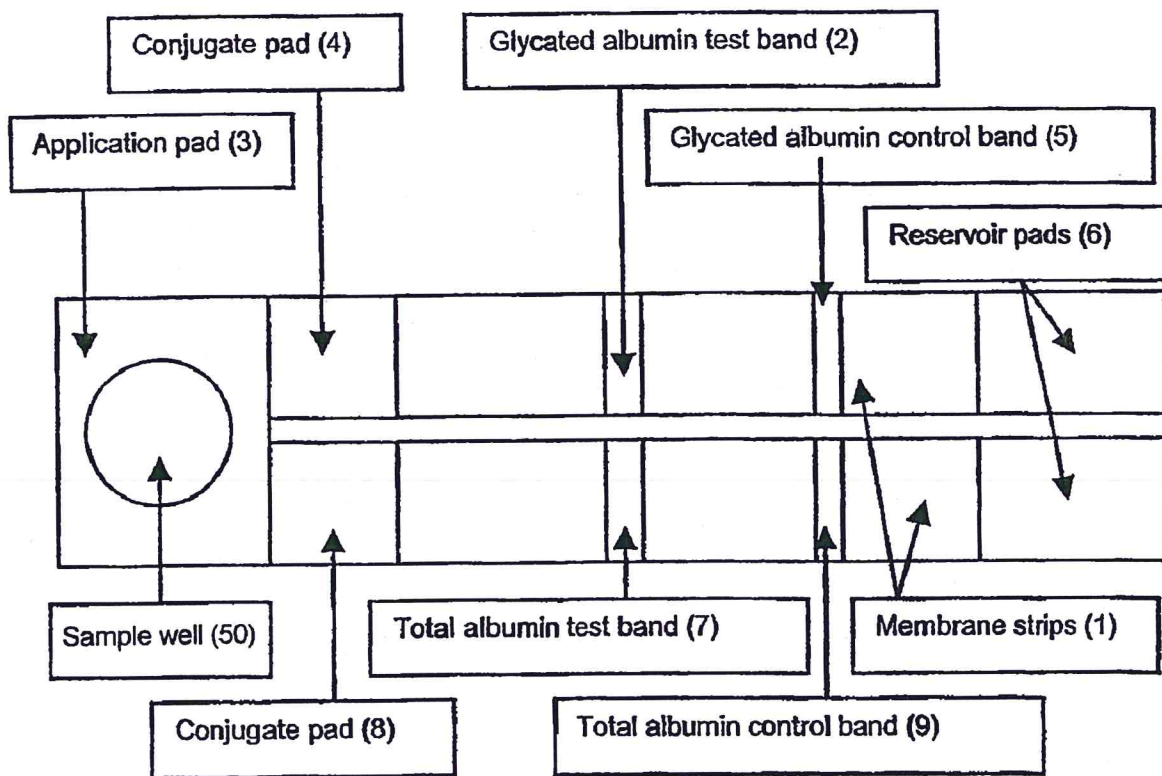

Test Strips: Overhead View

(57) **Abrégé/Abstract:**

A rapid immunochromatographic assay system is provided for measuring both glycated albumin and total albumin. The immunochromatographic assay system is comprised of a cassette that contains the test strips and testing reagents and a



(57) Abrégé(suite)/Abstract(continued):

measurement device that automatically reads, calculates and displays the test results. The test strip comprises a solid phase support (1) to which antibody to glycosylated albumin has been fixed to the support as a band (2). A sample application pad (3) contacts a conjugate pad (4) containing microparticles coated with anti-glycosylated albumin antibody. A control band (5) is provided to bind excess unreacted microparticles. A reservoir pad (6) is provided at the distal end to absorb excess sample fluid. The measurement device can be either a reflectance spectrometer or fluorometer. The level of glycosylated albumin is expressed as a percentage of the total albumin present in blood. The results of tests that are performed over a period of time are stored in the instrument's memory and presented in a numerical or graphical format so that an individual patient's glycosylated albumin level can be monitored over time.